

# XLS+ Multichannels



## **Multichannel XLS+ Pipettes**

---

Lighter weight

---

Less force

---

Greater control

---

Better balance

## **Lightweight Multichannel Pipettes**

Performance You Can Feel

**METTLER TOLEDO**

# Rainin XLS+ Multichannels

## Lightweight, Superior Performance

**The new Rainin XLS+™ multichannel pipettes by METTLER TOLEDO raise the bar on multichannel pipetting and performance. The new models are up to 35% lighter, which means they are easier to maneuver and easy on your hands, even after hours of pipetting.**

XLS+ pipettes are Rainin's most technically and ergonomically advanced multichannels available. Equipped with Rainin's patented LiteTouch™ System (LTS™) technology, these multichannel pipettes are easy on your hands and deliver precise and consistent sample pickup across all channels. Built rugged and with premium quality materials, these multichannels will deliver years of consistent, reliable service. And like all Rainin XLS+ pipettes, the embedded RFID tag supports advanced calibration tracking.




### Lightweight

At last... comfort and performance, even after hours of pipetting! Rainin XLS+ multichannels are lightweight, yet balanced and substantial, with a fit, feel and finish that is truly world class.



### Greater control

Rainin XLS+ multichannels are among the easiest and the most responsive multichannel pipettes you'll ever use. On all models and volume ranges, operation is smooth and absolutely consistent from start to finish.

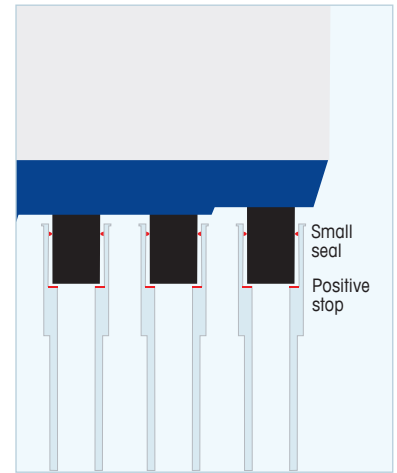

### Lightest touch

Choose Rainin's patented LiteTouch System (LTS) for ergonomics and ease of use. With a small seal and positive stop, tips load quickly, seal perfectly and are easy to eject.



**Channel consistency**

With a small seal and positive stop, LTS provides a fast, perfect seal across all channels. No rocking or hand-tightening required... ever! Perfect seals result in consistent sample pickup on all channels – first time, every time.



**Faster throughput**

Speed plate work with the Rainin E4 XLS+ electronic multichannel pipettes. Specialized modes streamline how you transfer multiple aliquots, perform serial dilutions and program complex pipetting routines.



**Pipet-Lite XLS+ and E4 XLS+ models**



# Rainin XLS+ Multichannel Pipettes

## Manual and Electronic Models



### Pipet-Lite XLS+ Multichannel

MT-Order No.	Description	Volume Range
<b>8-Channel</b>		
17013802	L8-10XLS+	0,5–10 µL
17013803	L8-20XLS+	2–20 µL
17013804	L8-50XLS+	5–50 µL
17013805	L8-200XLS+	20–200 µL
17013806	L8-300XLS+	20–300 µL
17014496	L8-1200XLS+	100–1200 µL
<b>12-Channel</b>		
17013807	L12-10XLS+	0,5–10 µL
17013808	L12-20XLS+	2–20 µL
17013809	L12-50XLS+	5–50 µL
17013810	L12-200XLS+	20–200 µL
17013811	L12-300XLS+	20–300 µL
17014497	L12-1200XLS+	100–1200 µL

### E4 XLS+ Multichannel

MT-Order No.	Description	Volume Range
<b>8-Channel</b>		
17013792	E8-10XLS+	0,5–10 µL
17013793	E8-20XLS+	2–20 µL
17013794	E8-50XLS+	5–50 µL
17013795	E8-200XLS+	20–200 µL
17013796	E8-300XLS+	20–300 µL
17014498	E8-1200XLS+	100–1200 µL
<b>12-Channel</b>		
17013797	E12-10XLS+	0,5–10 µL
17013798	E12-20XLS+	2–20 µL
17013799	E12-50XLS+	5–50 µL
17013800	E12-200XLS+	20–200 µL
17013801	E12-300XLS+	20–300 µL
17014499	E12-1200XLS+	100–1200 µL

Rainin LTS pipettes are designed for use with LTS tips.

RAININ TRADEMARKS: Rainin, Pipet-Lite, E4, XLS, XLS+, LTS, LiteTouch and Pipetting 360° are trademarks of Mettler-Toledo Rainin, LLC.

[www.mt.com/rainin](http://www.mt.com/rainin)

For more information

#### Mettler-Toledo Rainin, LLC

7500 Edgewater Drive, Oakland, CA 94621

Phone +1 510 564 1600

Fax +1 510 564 1604

Subject to technical changes

© 03/2016 Mettler-Toledo Rainin, LLC

Printed on demand from online content. 17700754 Rev E